

COSPIXZEL-1 Android System

Fine Lines Removal ▪ Acne Scars Removal ▪ Deep Wrinkles ▪ Face Tightening ▪ Stretch Marks
Fine Lines Removal ▪ Acne Scars Removal ▪ Deep Wrinkles ▪ Face Tightening ▪ Stretch Marks



Fine Lines
Removal



Acne Scars
Removal




Face
Tightening



Stretch
Marks



Deep
Wrinkles

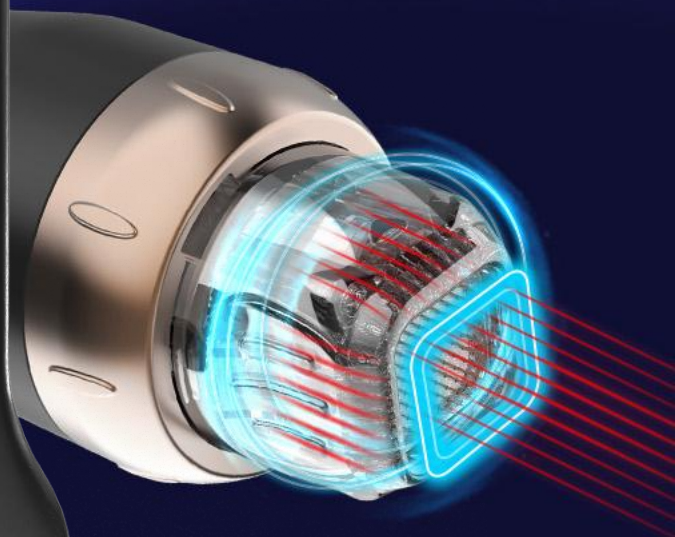



Application

One device. All skin types. All skin resurfacing

Face

Fine lines



Deep Wrinkles



Acne scars



Tightening



Skin reconstruction



Appearance of large pores



Body

Stretch Marks



Scars from injury or surgery



Fine lines






Advantage - With Vacuum

- MICRO NEEDLE RF (RF) applicator for procedures requiring ablation and resurfacing of the skin. Add vacuum - suction combined probe for better skin contact, especially around eyes and neck areas.




M-18

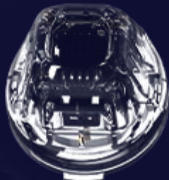

M-25



M-49

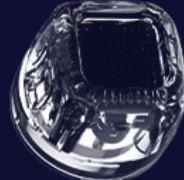

M-81

Stepping Motor Type

- Stepping Motor Type, different to existing solenoid type, the needle inserts to the skin smoothly with no shock, and causes no bleeding and no pain after the procedure

Needle Depth Control



- 0.4-4.0mm [0.2mm Step]
- Operate epidermis layer and dermis layer by controlling the needle depth in unit of 0.2mm.

Gold Plating

- Needle is durable and also have high Biocompatibility by applying Gold Plating. Patient with Metal allergy could also use it with not concerning Contact Dermatitis.


Seafy Needle System/Needle Thickness: 0.2mm

- Sterilized disposable needle tip
- operator can easily notice the applying RF energy from LED light.
- The needle structure is easy to penetrate the skin with minimum resistance.

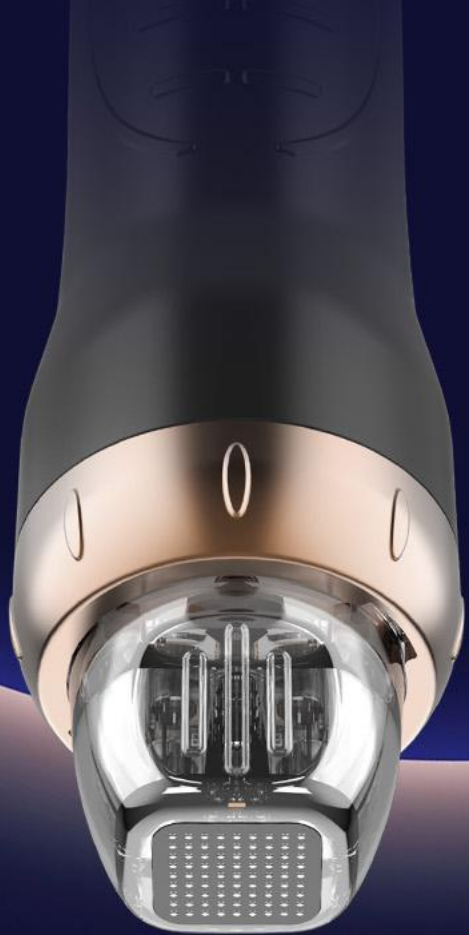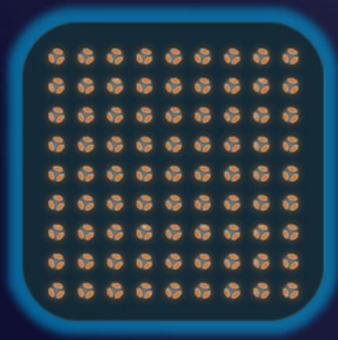

COSPIXZEL-1 -Fractional RF


FRACTIONAL RF applicator for non-invasive wrinkle reduction ,
suit for the patients fear of needle, recovery soon

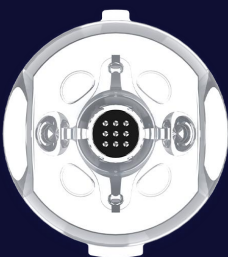

F-3*3


■ Eyes

■ Nose Area

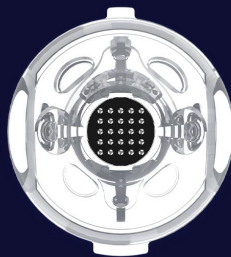

F-5*5

■ Face Area

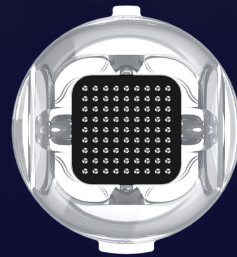

F-9*9

■ Body Area

COSPIXZEL-1-Microneedle RF


-- Negative polar (MONO/BI POLAR)

Negative Handle:


Together with MONO polar,
energy transfer Deeper

Under Bipolar

50% percent needle is positive and the other 50% is negative polar. this make the energy only works on depth of the needle. This mode is more suitable for face treatment.

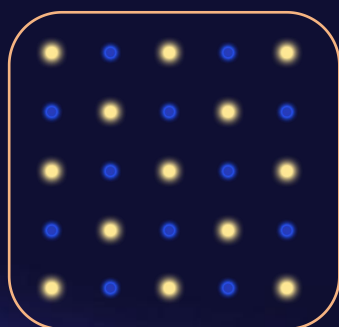

Face : Evenly

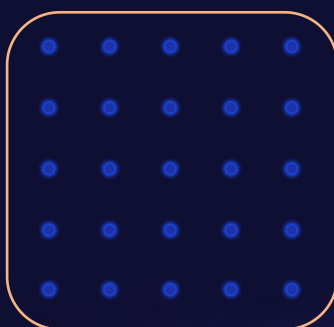

Body: Deeper

Under the Mono mode

all the needles on the tip is positive polar. The customer need to hold the negative polar in hand during the treatment. The RF energy can go to more deeper to skin than Bipolar mode. So it's more suitable for body treatment.

Cooling handle:

Max Reduce Pain


■ NEGATIVE POLE ■ POSITIVE ELECTRODE



- Beauticians administration
- Patients administration

Clients information storage

Treatment parameter date memory system

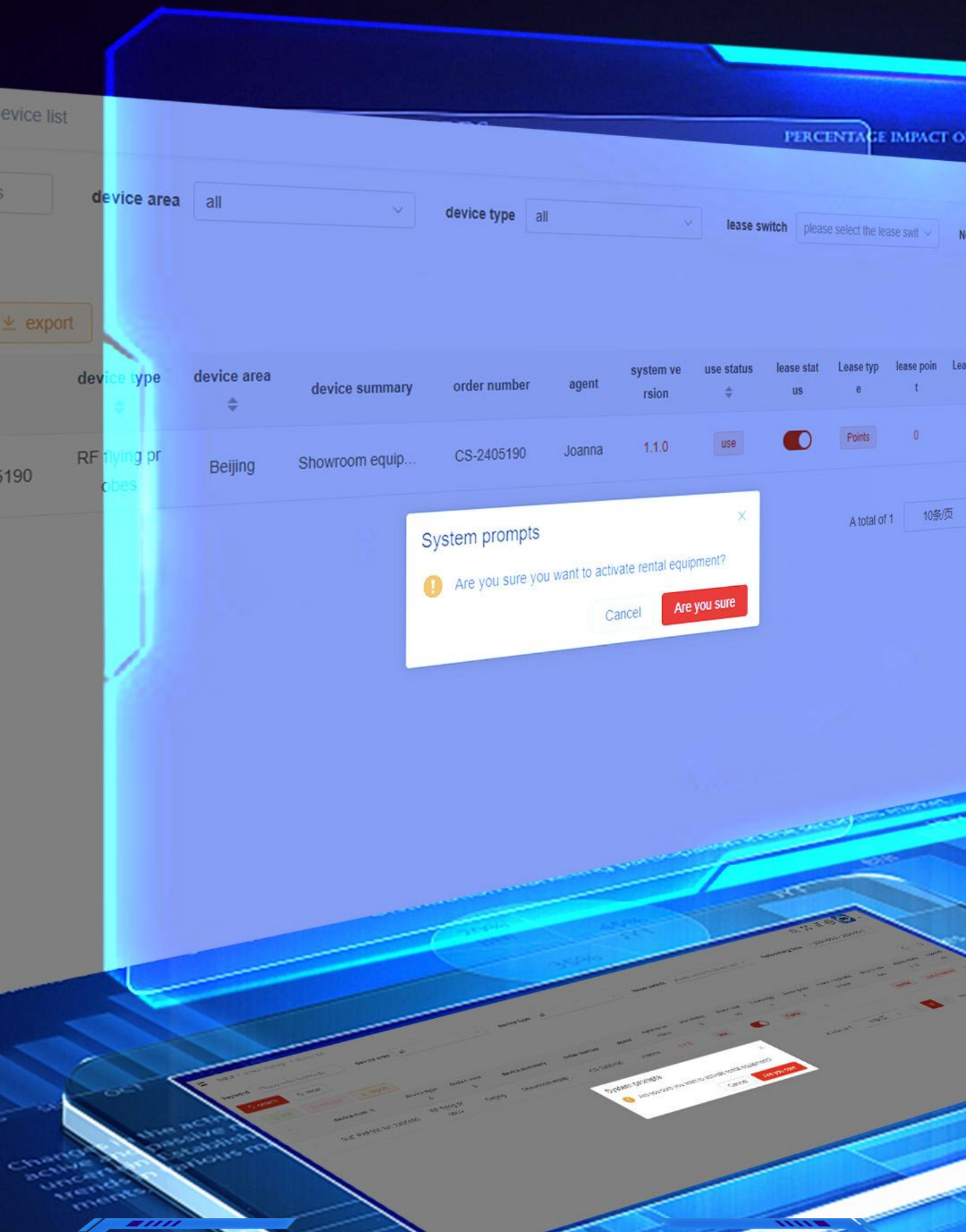

APP or PC control Rent-System

Distributers control the machine Anytime Anywhere

Shots limited or period time renting optional

Android smart system

WIFI

Bluetooth


Bluetooth

with Bluetooth,operator can control the machine
with computer or Mobile Phone.

WIFI


Machine connect with WIFI,record machine working
status and upload to Amazon Web Service.Distributor
remote control machine

vacuum turn on

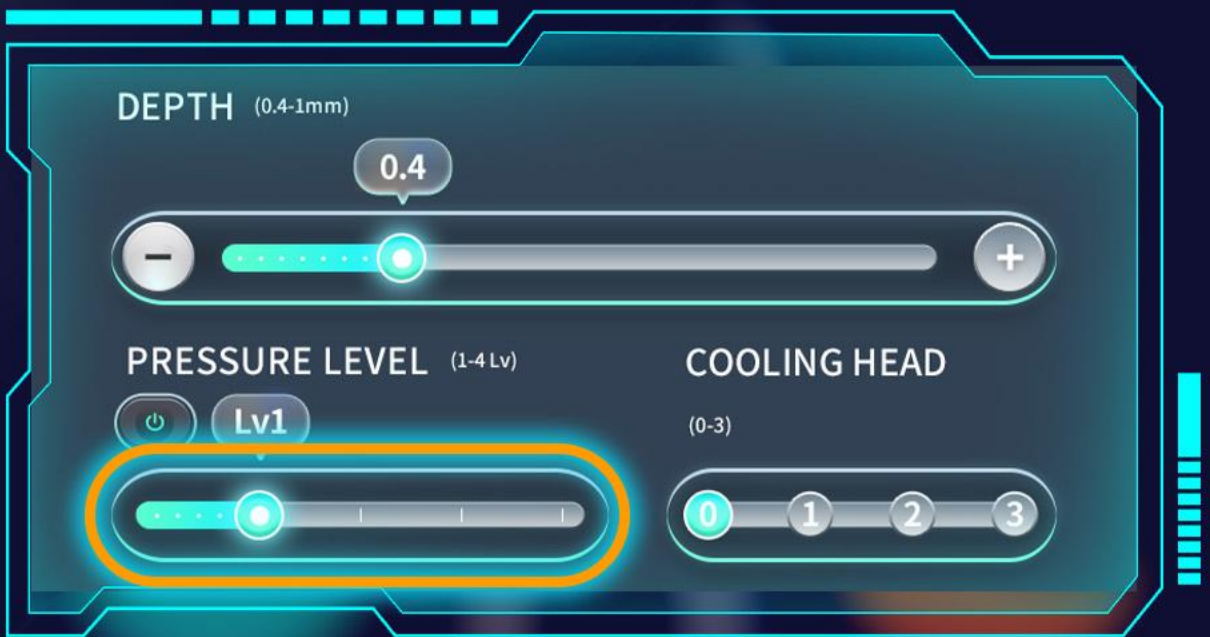

When the negative pressure is on, the needle length is 0.4-4mm adjustable

vacuum turn off

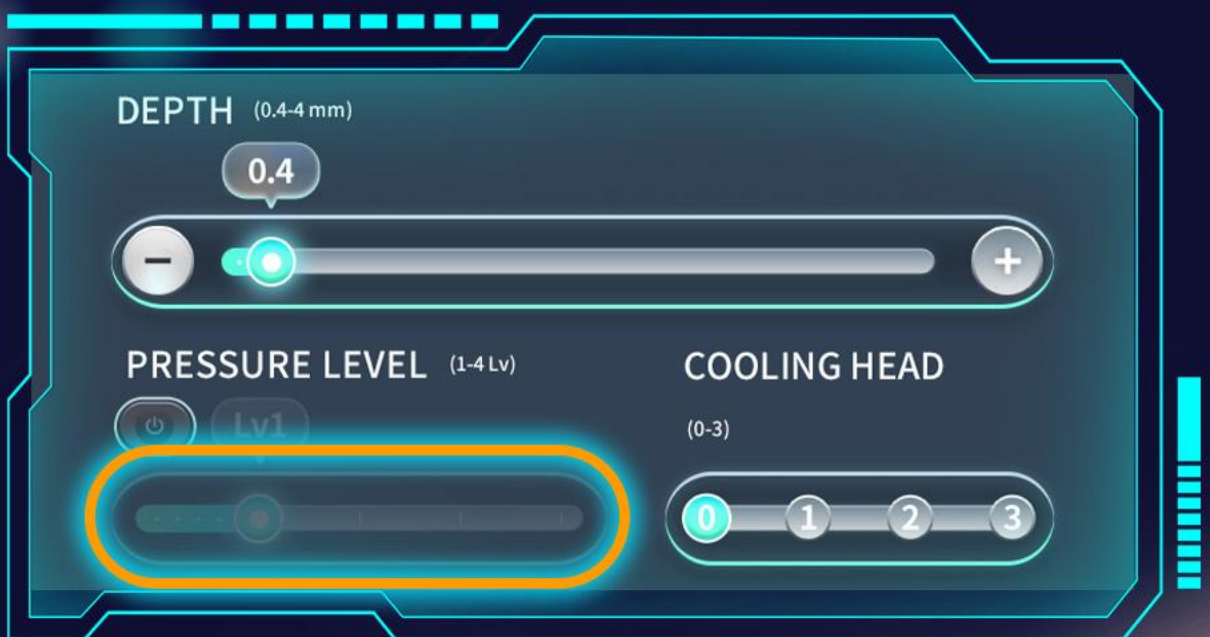

When closed, the needle length is adjustable from 0.4 to 1

The background features a series of concentric dotted circles in shades of teal and blue, creating a tunnel-like effect. In the center, there is a faint grid pattern. The word "THANKS" is prominently displayed in the middle of the grid.

THANKS